abcam

Product datasheet

Anti-E2F1 (phospho T433) antibody ab55325

★★★★★ 3 Abreviews 5 References 2 Images

Overview

Product name Anti-E2F1 (phospho T433) antibody

Description Rabbit polyclonal to E2F1 (phospho T433)

Host species Rabbit

Tested applications Suitable for: WB, IHC-P

Species reactivity Reacts with: Human

Predicted to work with: Mouse

Immunogen Synthetic peptide corresponding to Human E2F1 aa 350-450 (phospho T433).

Database link: Q01094

Positive control HeLa cell extract.

General notesThe Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or

contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

Properties

Form Liquid

Storage instructions Shipped at 4°C. Store at -20°C. Stable for 12 months at -20°C.

Storage buffer pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: 0.87% Sodium chloride, 50% Glycerol (glycerin, glycerine), PBS

Purity Immunogen affinity purified

Purification notesThe antibody against non-phosphopeptide was removed by chromatography using non-

phosphopeptide corresponding to the phosphorylation site.

Clonality Polyclonal

Isotype IgG

1

Applications

The Abpromise guarantee

Our <u>Abpromise guarantee</u> covers the use of ab55325 in the following tested applications.


The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

Application	Abreviews	Notes
WB		1/500 - 1/1000. Detects a band of approximately 47 kDa (predicted molecular weight: 47 kDa).
IHC-P	★★★★ (1)	Use at an assay dependent concentration.

Target

Function	Transcription activator that binds DNA cooperatively with dp proteins through the E2 recognition site, 5'-TTTC[CG]CGC-3' found in the promoter region of a number of genes whose products are involved in cell cycle regulation or in DNA replication. The DRTF1/E2F complex functions in the control of cell-cycle progression from G1 to S phase. E2F-1 binds preferentially RB1 protein, in a cell-cycle dependent manner. It can mediate both cell proliferation and p53-dependent apoptosis.	
Sequence similarities	Belongs to the E2F/DP family.	
Post-translational modifications	Phosphorylated by CDK2 and cyclin A-CDK2 in the S-phase. Acetylation stimulates DNA-binding. Enhanced under stress conditions such as DNA damage and inhibited by retinoblastoma protein pRB. Regulated by KAP1/TRIM28 which recruits HDAC1 to E2F1 resulting in deacetylation. Acetylated by P/CAF/KAT2B.	
Cellular localization	Nucleus.	

Images

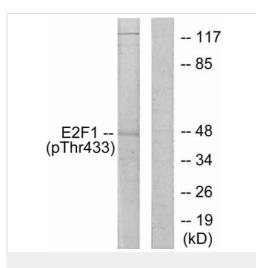

Western blot - Anti-E2F1 (phospho T433) antibody (ab55325)

All lanes : Anti-E2F1 (phospho T433) antibody (ab55325) at 1/500 dilution

Lane 1 : HeLa cell extract treated with Etoposide (at $25\mu M$ for 24 hrs).

Lane 2 : HeLa cell extract treated with Etoposide (at $25\mu M$ for 24 hrs), and with the immunising phosphopeptide.

Predicted band size: 47 kDa **Observed band size:** 47 kDa

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-E2F1 (phospho T433) antibody (ab55325)

Ab55325 staining human normal pancreas. Staining is localised to nuclear compartment.

Left panel: with primary antibody at 4 ug/ml. Right panel: isotype control.

Sections were stained using an automated system DAKO Autostainer Plus, at room temperature. Sections were rehydrated and antigen retrieved with the Dako 3-in-1 antigen retrieval buffer citrate pH6 in a DAKO PT Link. Slides were peroxidase blocked in 3% H2O2 in methanol for 10 minutes. They were then blocked with Dako Protein block for 10 minutes (containing casein 0.25% in PBS) then incubated with primary antibody for 20 minutes and detected with Dako Envision Flex amplification kit for 30 minutes. Colorimetric detection was completed with diaminobenzidine for 5 minutes. Slides were counterstained with Haematoxylin and coverslipped under DePeX. Please note that for manual staining we recommend to optimize the primary antibody concentration and incubation time (overnight incubation), and amplification may be required.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.com/abpromise or contact our technical team.

Terms and conditions

. Guarantee only valid for products bought direct from Abcam or one of our authorized distributors