

Product datasheet

Annexin V-iFluor 700 Apoptosis Staining / Detection Reagent ab219910

1 Image

Overview

Product name	Annexin V-iFluor 700 Apoptosis Staining / Detection Reagent
Detection method	Fluorescent
Sample type	Adherent cells, Suspension cells
Assay type	Quantitative
Species reactivity	Reacts with: Mammals, Other species
Product overview	Annexin V-iFluor 700 Apoptosis Staining / Detection Reagent (ab219910) is a cell- impermeable reagent designed to bind to phosphatidylserine (PS) residues exposed on the outer cell surface of cells with a flow cytometer or fluorescence microscopy at Ex/Em = 690/716 nm.

We recommend using an impermeable nuclear stain such as Propidium Iodide ([ab14083](#)) together with Annexin V-iFluor 700 Detection Reagent to discriminate necrotic and dead cells: plasma membrane is disrupted in these cells and therefore the Annexin V reagent will bind to PS found in the interior of cells.

Notes

Apoptosis is a regulated process of cell death that occurs during embryonic development as well as maintenance of tissue homeostasis. Inappropriately regulated apoptosis is implicated in different disease states, such as neurodegeneration disease and cancer. The apoptosis program is characterized by morphologic features, including loss of plasma membrane asymmetry and attachment, condensation off the cytoplasm and nucleus, and compaction and fragmentation of the nuclear chromatin. Exposure of phosphatidylserine (PS) on the external surface of the cell membrane has been reported to occur in the early phases of apoptotic cell death, during which the cell membrane remains intact. In leukocyte apoptosis, PS on the outer surface of the cell marks the cell for recognition and phagocytosis by macrophages. The human vascular anticoagulant, annexin V, is a 35-36 kDa Ca^{2+} dependent phospholipid binding protein that has a high affinity for PS, and shows minimal binding to phosphatidylcholine and sphingomyelin. Changes in PS asymmetry, which can be analyzed by measuring annexin V binding to the cell membrane, are generally observed before morphological changes associated with apoptosis occurred and before membrane integrity is lost.

Platform Flow cytometer

Properties

Storage instructions

Store at -20°C. Please refer to protocols.

Components	100 tests
Annexin V-iFluor 700 conjugate	1 x 200µl

Function

This protein is an anticoagulant protein that acts as an indirect inhibitor of the thromboplastin-specific complex, which is involved in the blood coagulation cascade.

Involvement in disease

Pregnancy loss, recurrent, 3

Sequence similarities

Belongs to the annexin family.
Contains 4 annexin repeats.

Domain


The [IL]-x-C-x-x-[DE] motif is a proposed target motif for cysteine S-nitrosylation mediated by the iNOS-S100A8/A9 transnitrosylase complex.

A pair of annexin repeats may form one binding site for calcium and phospholipid.

Post-translational modifications

S-nitrosylation is induced by interferon-gamma and oxidatively-modified low-density lipoprotein (LDL(ox)) possibly implicating the iNOS-S100A8/9 transnitrosylase complex.

Images

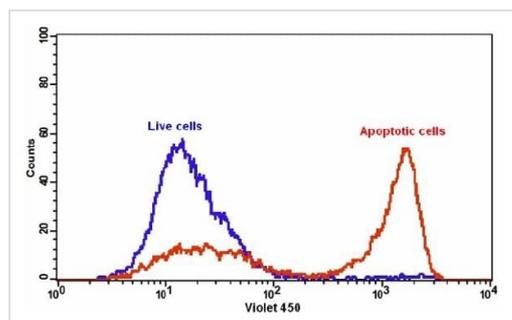

The detection of binding activity of

[ab219911](#)

Annexin V-mFluor Violet450 Detection Reagent ([ab219911](#)), which has similar PS binding properties to Annexin V-iFluor 700 conjugate. Detection of phosphatidylserine (PS) exposure in Jurkat cells. Jurkat cells were left untreated (blue) or treated with 1 µM staurosporine (red) in a 37°C, 5% CO₂ incubator for 5 hours. Cells were then incubated with Annexin V Detection Reagent for 30 minutes. The fluorescence intensity of Annexin V-mFluor 450 Detection Reagent was measured with a FACSCalibur (BD systems) flow cytometer using violet laser at Ex/Em = 405/450 nm.

In live non-apoptotic cells, Annexin V-iFluor 700 conjugate detects innate apoptosis in non-induced cells, which is typically 2-6% of all cells. In apoptotic cells Annexin V-iFluor 700 conjugate binds to phosphatidylserine, which is located on the outer leaflet of the cell membrane, resulted in increased staining intensity.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Promise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you

- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors