

Product datasheet

Eosin Y Solution (Modified Alcoholic) ab246824

1 Image

Overview

---

<b>Product name</b>	Eosin Y Solution (Modified Alcoholic)
<b>General notes</b>	<p>Eosin Y Solution (Modified Alcoholic) ab246824 is intended for use in histology and cytology applications. This newly formulated Eosin provides the benefits of a traditional alcoholic formula with significant improvements in usability. Advantages include lower evaporation rate, better color patterns, reduced tendency to spill over container, hands, and countertops, and improved surface tension to remain on tissue section.</p> <p>Cytoplasm: Light Pink            Collagen: Pink            Muscle: Pink/Rose            Erythrocytes: Pink/Red</p> <p><b>Control Tissue:</b> Any well fixed paraffin embedded or frozen tissue section. Cell smear.</p> <p><b>Procedure:</b></p> <ol style="list-style-type: none"> <li>1. Deparaffinize sections if necessary and hydrate to distilled water.</li> <li>2. Dip slide in absolute alcohol and blot excess off.</li> <li>3. Apply adequate Eosin Y Solution (Modified Alcoholic) to completely cover tissue section to excess and incubate for 2-3 minutes.</li> <li>4. Rinse slide using absolute alcohol.</li> <li>5. Dehydrate slide in 3 changes of absolute alcohol.</li> <li>6. Clear slide and mount in synthetic resin.</li> </ol> <p><b>Other products for staining tissue sections</b></p> <p>This Eosin Y solution is also available in <a href="#">H&amp;E staining kit ab245880</a>.</p> <p>Find more kits and reagents in the <a href="#">special stains guide</a>, or products for antigen retrieval, blocking, signal amplification, visualization, counterstaining, and mounting in the <a href="#">IHC kits and reagents guide</a>.</p>


Properties

---

<b>Form</b>	Liquid
<b>Storage instructions</b>	Shipped at Room Temperature. Store at Room Temperature.
<b>Storage buffer</b>	Constituents: 65% Ethanol, 3.8% Isopropanol, 1% Acetic acid, 0.25% Eosin Y

## Images

---

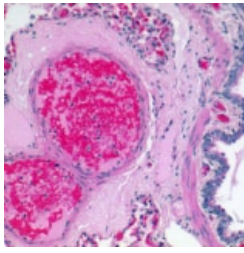

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Eosin Y Solution (Modified Alcoholic) (ab246824)

ab246824 Eosin Y Solution (Modified Alcoholic) staining formalin-fixed-paraffin embedded cell smear.

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

### Our Abpromise to you: Quality guaranteed and expert technical support

---

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

### Terms and conditions

---

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors