

Product datasheet

Anti-ICOS antibody [SP98] ab105227

Recombinant RabMAb

★ ★ ★ ☆ ☆ 2 Abreviews 10 References 8 Images

Overview

Product name	Anti-ICOS antibody [SP98]
Description	Rabbit monoclonal [SP98] to ICOS
Host species	Rabbit
Tested applications	Suitable for: IHC-P, IHC-Fr
Species reactivity	Reacts with: Human Predicted to work with: Mouse
Immunogen	Synthetic peptide within Human ICOS aa 150 to the C-terminus (C terminal). The exact sequence is proprietary. Database link: Q9Y6W8
Positive control	IHC-P: Human lymph node, pancreatic adenocarcinoma, thymus and tonsil tissue. IHC-Fr: Frozen human tonsil tissue sections
General notes	This product is a recombinant monoclonal antibody, which offers several advantages including: <ul style="list-style-type: none"> - High batch-to-batch consistency and reproducibility - Improved sensitivity and specificity - Long-term security of supply - Animal-free production For more information see here . This product is FOR RESEARCH USE ONLY. For commercial use, please contact partnerships@abcam.com.

Properties

Form	Liquid
Storage instructions	Shipped at 4°C. Upon delivery aliquot and store at -20°C. Avoid freeze / thaw cycles.
Storage buffer	pH: 7.60 Preservative: 0.1% Sodium azide Constituents: PBS, 1% BSA
Purity	Protein A/G purified
Purification notes	Purified from TCS by protein A/G.

Clonality	Monoclonal
Clone number	SP98
Isotype	IgG

Applications


The Abpromise guarantee Our [Abpromise guarantee](#) covers the use of ab105227 in the following tested applications. The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

Application	Abreviews	Notes
IHC-P	★★★★★ (1)	1/100. Perform antigen retrieval before commencing with IHC staining protocol by boiling tissue section in EDTA buffer, pH 8.0 for 10 min followed by cooling at RT for 20 min.
IHC-Fr	★☆☆☆☆ (1)	Use a concentration of 5 µg/ml.

Target


Function	Enhances all basic T-cell responses to a foreign antigen, namely proliferation, secretion of lymphokines, up-regulation of molecules that mediate cell-cell interaction, and effective help for antibody secretion by B-cells. Essential both for efficient interaction between T and B-cells and for normal antibody responses to T-cell dependent antigens. Does not up-regulate the production of interleukin-2, but superinduces the synthesis of interleukin-10. Prevents the apoptosis of pre-activated T-cells. Plays a critical role in CD40-mediated class switching of immunoglobulin isotypes.
Tissue specificity	Activated T-cells. Highly expressed on tonsillar T-cells, which are closely associated with B-cells in the apical light zone of germinal centers, the site of terminal B-cell maturation. Expressed at lower levels in thymus, lung, lymph node and peripheral blood leukocytes. Expressed in the medulla of fetal and newborn thymus.
Involvement in disease	Immunodeficiency, common variable, 1
Sequence similarities	Contains 1 Ig-like V-type (immunoglobulin-like) domain.
Post-translational modifications	N-glycosylated.
Cellular localization	Secreted and Cell membrane.

Images

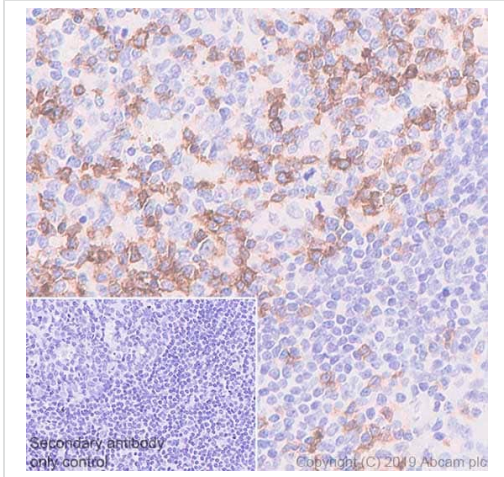


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-ICOS antibody [SP98] (ab105227)

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) analysis of Human tonsil tissue sections labeling ICOS with ab105227 at 1/100 dilution (0.28 µg/ml). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0, epitope retrieval solution 2) for 30 mins. Rabbit specific IHC polymer detection kit HRP/DAB (ab209101) was used as the secondary antibody. Hematoxylin was used as a counterstain. Sporadically positive staining in the germinal center of human tonsil, performed on a Leica Biosystems BOND™ RX instrument. The section was incubated with ab105227 for 30 mins at room temperature.

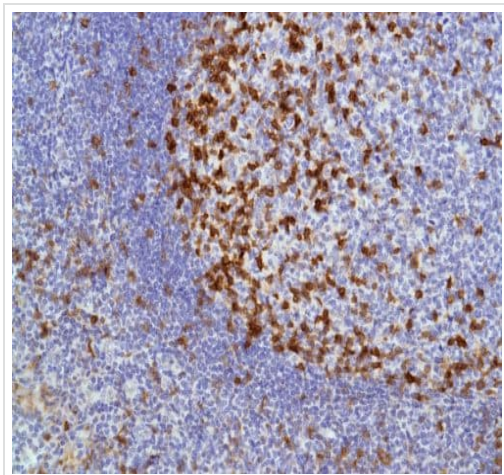


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-ICOS antibody [SP98] (ab105227)

Formalin-fixed, paraffin-embedded human tonsil tissue stained for ICOS using ab105227 at 1/100 dilution in immunohistochemical analysis.

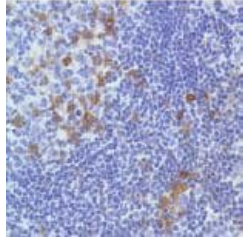


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-ICOS antibody [SP98] (ab105227)

ab105227, at 1/100 dilution, staining ICOS in formalin-fixed, paraffin-embedded Human reactive lymph node tissue by Immunohistochemistry.

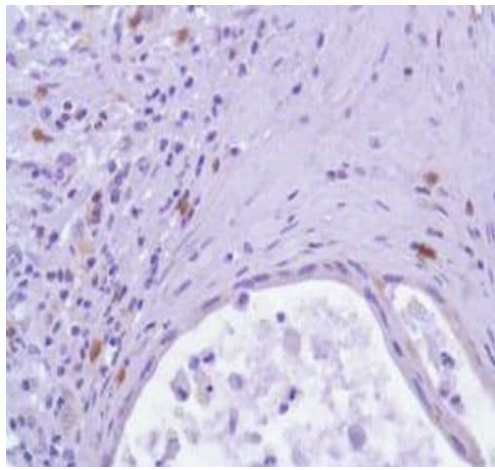


Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-ICOS antibody [SP98] (ab105227)

Formalin-fixed, paraffin-embedded human pancreatic adenocarcinoma tissue stained for ICOS using ab105227 at 1/100 dilution in immunohistochemical analysis.

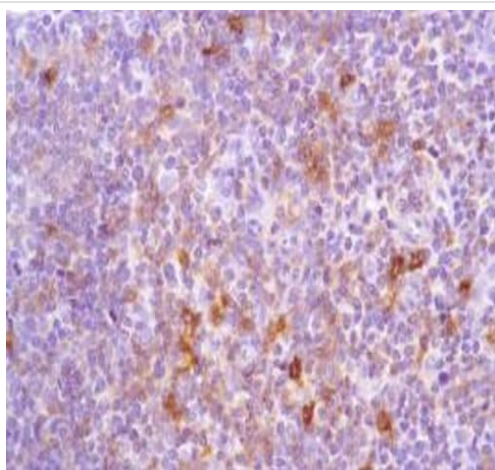

Immunohistochemistry (Formalin/PFA-fixed paraffin-embedded sections) - Anti-ICOS antibody [SP98] (ab105227)

Formalin-fixed, paraffin-embedded human thymus tissue stained for ICOS using ab105227 at 1/100 dilution in immunohistochemical analysis.

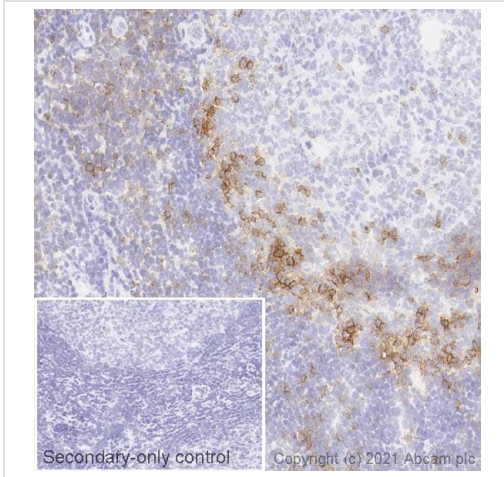


Immunohistochemistry (Frozen sections) - Anti-ICOS antibody [SP98] (ab105227)

IHC image of ICOS staining in a section of frozen normal human tonsil* performed on a Leica BOND™ system using the standard protocol. The section was fixed in 10% paraformaldehyde (10 min) prior to staining. The section was incubated with ab105227, 5ug/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset secondary-only control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

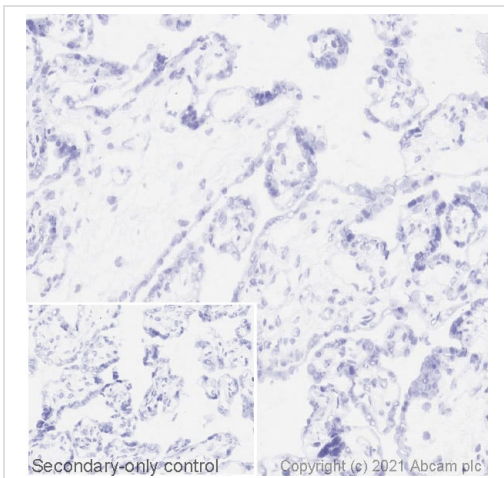

Immunohistochemistry (Frozen sections) - Anti-ICOS antibody [SP98] (ab105227)

Negative control image: IHC image of ICOS staining in a section of frozen normal human placenta* performed on a Leica BOND™ system using the standard protocol. The section was fixed in 10% paraformaldehyde (10 min) prior to staining. The section was incubated with ab105227, 5ug/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. The inset secondary-only control image is taken from an identical assay without primary antibody.

For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times.

*Tissue obtained from the Human Research Tissue Bank, supported by the NIHR Cambridge Biomedical Research Centre

Why choose a recombinant antibody?

Research with confidence
Consistent and reproducible results

Long-term and scalable supply
Recombinant technology

Success from the first experiment
Confirmed specificity

Ethical standards compliant
Animal-free production

Anti-ICOS antibody [SP98] (ab105227)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors