

Product datasheet

InstantBlue® Coomassie Protein Stain (ISB1L) ab119211

★★★★★ 4 Abreviews 492 References 2 Images

Overview

Product name InstantBlue® Coomassie Protein Stain (ISB1L)

General notes InstantBlue® is a ready to use Coomassie protein stain for polyacrylamide gels. Its unique mechanism of action stains proteins in 15 minutes, while leaving a clear background eliminating the need to fix, wash or destain. Formulated for safe use and easy disposal, it's ready to use straight out of the bottle and comes in convenient premixed one liter bottles.

Why InstantBlue® Protein Stain?

There is no need to wash, fix, microwave or destain gels when using InstantBlue™. Homemade stains and premade solutions require several or all of these extra steps. InstantBlue uses a simple one step protocol.

The unique patented mechanism for rapid Coomassie blue staining of protein gels begins in moments, and results are achieved within 15 minutes. Some premade and traditional homemade Coomassie R-250 protein stains can take three hours or more to fully stain gels, and then require destaining typically overnight.

Unparalleled sensitivity allows protein levels as low as 5ng per band to be detected (BSA) when staining overnight. This is thanks to InstantBlue's protein specific staining, which gives a clear background, improving the signal / noise ratio of the stain. Despite the staining power, gels can be left in InstantBlue for weeks and still come out with a crystal clear background, all without the need for destaining or protein stain remover.

The nontoxic formula is sink disposable and thanks to its methanol free composition will not shrink gels. Other stains require solvent disposal procedures and can require microwaving, which produces dangerous fumes. The effective staining and safe composition of InstantBlue means that there's no need for microwaving, solvent disposal, or fume hoods.

InstantBlue will stain at least 40 mini gels per bottle, boasting economical staining with only 25ml per mini gel. As a result, this allows you to use less protein staining solution when compared to 50–100ml of stain per gel for homemade solutions and some premixed solutions.

Fully compatible with mass spectrometry, gels can be destained using standard procedures. The protein staining solution's acetic acid free formulation means it will not methylate or acetylate protein during the protein staining protocol.

Key Benefits

Ultra fast staining – Results in 15 minutes or less

One step procedure – No washing, fixing, microwaving or destaining

High sensitivity – 5ng bands detectable

Clear background – High signal / noise ratio

Flexible – No overstaining

Reliable – Stringent batch to batch consistency

Efficient – Only 25ml per gel

Safe composition – Nontoxic; no fume hood or solvent disposal required

Methanol free – No gel shrinkage or protein methylation

Acetic acid free – No protein acetylation

Mass spectrometry compatible – Destainable. No residual methylation or acetylation

This product is manufactured by Expedeon, an Abcam company. Expedeon product code ISB1L was previously called InstantBlue Protein Stain 1 L and is the same as the 1000 ml size of this product.


Manufactured under patent numbers : US 12/226,797 ; SG 201103010-3 ; SG 200807953-5 ; JP 2009-507182 ; IN 5722/CHENP/2008 ; GB 0608377.8 ; EP 10187154.9; EP 07733647.7 ; CN 200780021434.8 ; CA 2,649,835 ; AU 2007245443

INSTANTBLUE is a registered trademark of Abcam in the US and an unregistered trademark elsewhere


Properties

Form	Liquid
Storage instructions	Shipped at Room Temperature. Store at +4°C.
Storage buffer	Constituents: 0.56% Ethanol, 15.19% Phosphoric acid

Images

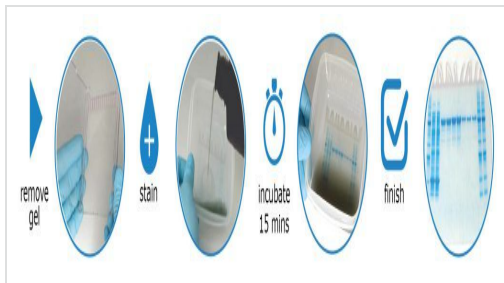

Western blot - InstantBlue® Coomassie Protein
Stain (ISB1L) (ab119211)

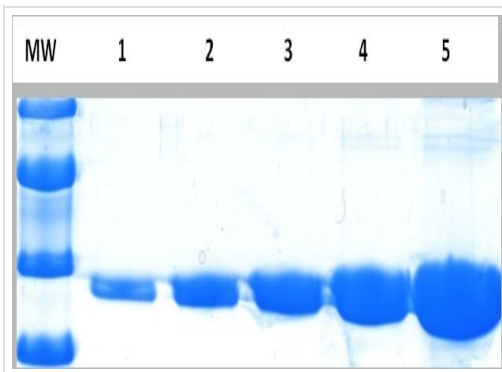

Western blot - InstantBlue® Coomassie Protein

Stain (ISB1L) (ab119211)

This image is courtesy of an anonymous abreview.

SDS-PAGE image showing molecular weight marker (MW) and His-tag recombinant protein fractions purified in Ni-NTA resin in lanes 1-5 stained with InstantBlue. InstantBlue was easy to handle and use - just pour a small amount onto your gel and within 10 minutes you can see the protein bands! No methanol and glacial acetic acid - no smell!

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors