

Product datasheet

Metformin hydrochloride, AMPK activator ab120847

4 References 2 Images


Overview

Product name	Metformin hydrochloride, AMPK activator
Description	AMPK activator
Biological description	AMPK activator. Antidiabetic and anticancer activity <i>in vivo</i> . Rapidly crosses the blood-brain barrier. Also available in simple stock solutions (ab146725) - add 1 ml of water to get an exact, ready-to-use concentration.

Purity > 99%

CAS Number 1115-70-4

Chemical structure

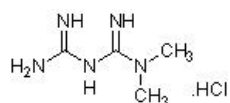

Properties


Chemical name	<i>N,N</i> -Dimethylimidodicarbonimidic diamide hydrochloride
Molecular weight	165.63
Molecular formula	C ₄ H ₁₁ N ₅ .HCl
PubChem identifier	14219
Storage instructions	Store at Room Temperature. The product can be stored for up to 12 months.
Solubility overview	Soluble in water to 100 mM and in DMSO to 50 mM
Handling	Wherever possible, you should prepare and use solutions on the same day. However, if you need to make up stock solutions in advance, we recommend that you store the solution as aliquots in tightly sealed vials at -20°C. Generally, these will be useable for up to one month. Before use, and prior to opening the vial we recommend that you allow your product to equilibrate to room temperature for at least 1 hour. Need more advice on solubility, usage and handling? Please visit our frequently asked questions (FAQ) page for more details.
SMILES	CN(C)C(=N)N=C(N)N.Cl
Source	Synthetic

Applications

The Abpromise guarantee Our [Abpromise guarantee](#) covers the use of ab120847 in the following tested applications. The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.


Application	Abreviews	Notes
Functional Studies		Use at an assay dependent concentration.

Images

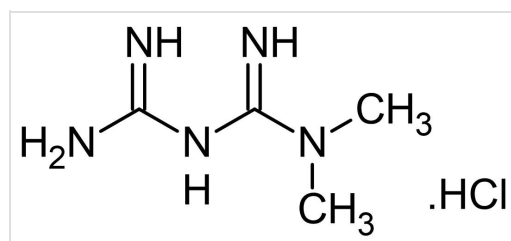

Chemical Structure - Metformin hydrochloride, AMPK activator (ab120847)

2D chemical structure image of ab120847, Metformin hydrochloride, AMPK activator

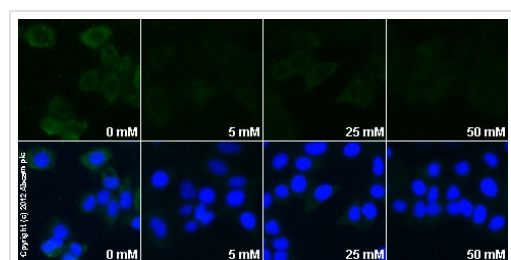

Functional Studies - Metformin hydrochloride, AMPK activator (ab120847)

[ab1](#) staining HIF1 α in MCF7 cells treated with metformin hydrochloride (ab120847), by ICC/IF. Decrease in HIF1 α expression correlates with increased concentration of metformin hydrochloride, as described in literature.

The cells were incubated at 37°C for 24h in media containing different concentrations of ab120847 (metformin hydrochloride) in water, fixed with 4% formaldehyde for 10 minutes at room temperature and blocked with PBS containing 10% goat serum, 0.3 M glycine, 1% BSA and 0.1% tween for 2h at room temperature. Staining of the treated cells with [ab1](#) (10 μ g/ml) was performed overnight at 4°C in PBS containing 1% BSA and 0.1% tween. A DyLight 488 goat anti-mouse polyclonal antibody ([ab96879](#)) at 1/250 dilution was used as the secondary antibody. Nuclei were counterstained with DAPI and are shown in blue.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES, NOT FOR USE IN HUMANS"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery

- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors
- Abcam biochemicals are novel compounds and we have not tested their biological activity in house. Please use the literature to identify how to use these products effectively. If you require further assistance please contact the scientific support team