

Human Cytokine Antibody Array F (40 Targets) - Quantitative ab197443

[3 Images](#)

Overview

Product name	Human Cytokine Antibody Array F (40 Targets) - Quantitative
Sample type	Serum, Plasma, Other biological fluids, Cell Lysate, Cell culture media, Tissue Lysate
Assay type	Quantitative
Assay time	6h 00m
Species reactivity	Reacts with: Human
Product overview	ab197443 can be used for quantitative measurement of 40 Human cytokines. Suitable for serum, plasma, cell culture media, other body fluids, cell and tissue lysates.

Targets names: ADAMTS13, Aggrecan, Angiotensinogen, B7-H1, BMPR-IA, BMPR-II, Cadherin-11, CD27, CD6, Ck beta 8-1, CNTF, DNAM-1, EMMPRIN, FLRG, Follistatin-like 1, Fractalkine, Galectin-1, GITR Ligand, Granulysin, IL-1 R3, IL-15 R, IL-17E, IL-32 alpha, L1CAM-2, LRIG3, LRP-6, MEPE, Nectin-4, Periostin, Persephin, Renin, RGM-B, ROBO3, S100A8, Siglec-7, Syndecan-3, Thrombospondin-2, Thrombospondin-5, Tie-1, ULBP-2.

[Learn more about quantitative antibody arrays](#)

Notes Quantitative antibody arrays can be used to quantitate up to 40 cytokines with as little as 50 µL of sample. Arrays are available for 400 human and 200 mouse proteins.

Each glass slide is spotted with 16 identical antibody arrays and is provided with a 16 well gasket to allow separate samples to be applied to each array. Each antibody is spotted in quadruplicate in each array. For reproducible quantitation, eight of the arrays are used with a cocktail of protein standards to produce a standard curve. The same 8 standard curve arrays can be used across multiple slides, allowing measurement of 22 experimental samples with 2 slides, 36 with 3 slides etc.

For high through-put, 4 slides can be nested into a tray matching a standard microplate, allowing automated processing with a liquid handling workstation.

Array processing can be completed within one working day (see workflow diagram below) and arrays can be analyzed with a wide number of laser glass slide array / gene microarray scanners. If you don't have a suitable scanner then we recommend our membrane antibody arrays, which

can be analyzed with any WB chemiluminescent reader.

[Learn more about quantitative antibody arrays](#)

[Learn more about membrane antibody arrays](#)

The standard is affinity purified.


Properties

Storage instructions


Store at -20°C. Please refer to protocols.

Components	8 tests	22 tests	50 tests
20X Wash Buffer I	2 x 30ml	3 x 30ml	6 x 30ml
20X Wash Buffer II	1 x 30ml	1 x 30ml	2 x 30ml
Cy3 equivalent dye-conjugated Streptavidin	1 x 5µl	2 x 5µl	4 x 5µl
Human Cytokine Array F (40 T) Biotinylated Antibody Cocktail	1 unit	2 units	4 units
Human Cytokine Array F (40 T) Glass Slide	1 slide	2 slides	4 slides
Human Cytokine Array F (40 T) Lyophilized Cytokine Standard Mix	1 vial	2 vials	2 vials
Sample Diluent	1 x 15ml	1 x 15ml	2 x 15ml
Slide Washer/Dryer (30 mL Centrifuge Tube)	1 unit	1 unit	2 units

Images

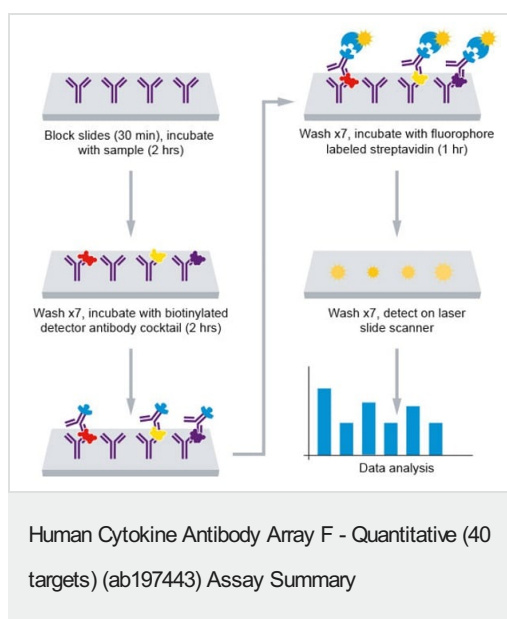

Assay summary for Abcam Human Cytokine Antibody Array F (40 Targets) - Quantitative (ab197443).

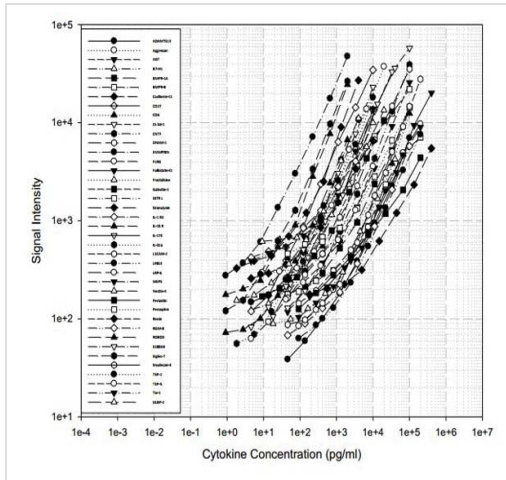

Standard Curve obtained using Abcam Human Cytokine Antibody Array F (40 Targets) - Quantitative (ab197443).

Human Cytokine Antibody Array F - Quantitative (40 targets) (ab197443) Standard Curve

POS1	POS2	ADAMTS13
Aggrecan	Angiotensinogen	B7-H1
BMPR-IA	BMPR-II	Cadherin-11
CD27	CD6	Ck beta 8-1
CNTF	DNAM-1	EMMPRIN
FLRG	Follistatin-like 1	Fractalkine
Galectin-1	GITR Ligand	Granulysin
IL-1 R3	IL-15 R	IL-17E
IL-32 alpha	L1CAM-2	LRIG3
LRP-6	MEPE	Nectin-4
Periostin	Persephin	Renin
RGM-B	ROBO3	S100A8
Siglec-7	Syndecan-3	Thrombospondin-2
Thrombospondin-5	Tie-1	ULBP-2

Array map for Abcam Human Cytokine Antibody Array F (40 Targets) - Quantitative (ab197443).

Human Cytokine Antibody Array F - Quantitative (40 targets) (ab197443) Array Map

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

Terms and conditions

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors