abcam

Product datasheet

Valproic acid, sodium salt, Histone deacetylase inhibitor ab120745

1 References 2 Images

Overview

Product name Valproic acid, sodium salt, Histone deacetylase inhibitor

Description Histone deacetylase inhibitor

Biological description Histone deacetylase inhibitor (IC $_{50}$ = 0.4 μ M for HDAC1). Increases GABA expression and

activates Wnt dependent gene expression. Centrally active upon oral and intravenous

administration.

CAS Number 1069-66-5

Chemical structure CO2Na

H₃C CH₃

Properties

Chemical name Sodium 2-propylpentanoate

Molecular weight 166.20

Molecular formula C₈H₁₅NaO₂

Storage instructions Store at Room Temperature. Store under desiccating conditions. The product can be stored for

up to 12 months.

Solubility overview Soluble in water to 100 mM

Handling Wherever possible, you should prepare and use solutions on the same day. However, if you need

to make up stock solutions in advance, we recommend that you store the solution as aliquots in tightly sealed vials at -20°C. Generally, these will be useable for up to one month. Before use, and

prior to opening the vial we recommend that you allow your product to equilibrate to room

temperature for at least 1 hour.

Toxic, refer to SDS for further information.

Need more advice on solubility, usage and handling? Please visit our frequently asked

questions (FAQ) page for more details.

Source Synthetic


The Abpromise guarantee

Our <u>Abpromise guarantee</u> covers the use of ab120745 in the following tested applications.


The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

| Application | Abreviews | Notes |
|--------------------|-----------|--|
| Functional Studies | | Use at an assay dependent concentration. |

Images

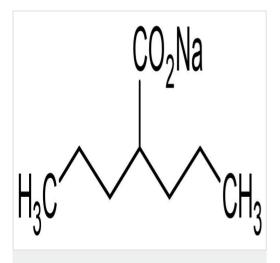

Chemical Structure - Valproic acid, sodium salt, Histone deacetylase inhibitor (ab120745) 2D chemical structure image of ab120745, Valproic acid, sodium salt, Histone deacetylase inhibitor

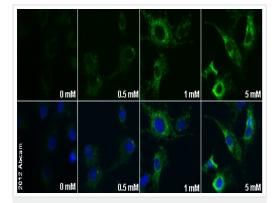

Immunocytochemistry/ Immunofluorescence -Valproic acid, sodium salt, Histone deacetylase inhibitor (ab120745) ab70362 staining adiponectin receptor 1 in HepG2 cells treated with valproic acid, sodium salt (ab120745), by ICC/IF. Increased cytoplasmatic staining of adiponectin receptor 1 correlates with increased concentration of valproic acid, as described in literature. The cells were incubated at 37°C for 24h in media containing different concentrations of ab120745 (valproic acid) in DMSO, fixed with 4% formaldehyde for 10 minutes at room temperature and blocked with PBS containing 10% goat serum, 0.3 M glycine, 1% BSA and 0.1% tween for 2h at room temperature. Staining of the treated cells with ab70362 (5 μg/ml) was performed overnight at 4°C in PBS containing 1% BSA and 0.1% tween. A DyLight 488 goat anti-rabbit polyclonal antibody (ab96899) at 1/250 dilution was used as the secondary antibody. Nuclei were counterstained with DAPI and are shown in blue.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES, NOT FOR USE IN HUMANS"

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.com/abpromise or contact our technical team.

Terms and conditions

- · Guarantee only valid for products bought direct from Abcam or one of our authorized distributors
- Abcam biochemicals are novel compounds and we have not tested their biological activity in house. Please use the literature to identify how to use these products effectively. If you require further assistance please contact the scientific support team