

# Cytokine Array - Mouse Cytokine Antibody Array (Membrane, 96 Targets) ab193659

[6 References](#) [1 Image](#)

### Overview

<b>Product name</b>	Cytokine Array - Mouse Cytokine Antibody Array (Membrane, 96 Targets)
<b>Sample type</b>	Serum, Plasma, Other biological fluids, Cell Lysate, Cell culture media, Tissue Lysate
<b>Assay type</b>	Semi-quantitative
<b>Species reactivity</b>	<b>Reacts with:</b> Mouse
<b>Product overview</b>	ab193659 can be used for the simultaneous detection of 96 Mouse Cytokine proteins in serum, plasma, cell culture media and other liquid samples types. It comes as 1 kit including 2 separate membranes labelled C3 and C4 (see below for images showing the array layout on each membrane).


**Target names:** Axl, Bfgf, BLC (CXCL13), CD30 Ligand (TNFSF8), CD30 (TNFRSF8), CD40 (TNFRSF5), CRG-2, CTACK (CCL27), CXCL16, CD26 (DPPIV), Dtk, Eotaxin-1 (CCL11), Eotaxin-2 (MPIF-2/CCL24), E-Selectin, Fas Ligand (TNFSF6), Fc gamma RIIb (CD32b), Flt-3 Ligand, Fractalkine (CX3CL1), GCSF, GTR (TNFRSF18), GM-CSF, HGFR, ICAM-1 (CD54) IFN-gamma, IGFBP-2, IGFBP-3, IGFBP-5, IGFBP-6, IGF-1, IGF-2, IL-1 beta (IL-1 F2), IL-10, IL-12 p40/p70, IL-12 p70, IL-13, IL-15, IL-17A, IL-17 RB, IL-1 alpha (IL-1 F1), IL-2, IL-3, IL-3 R beta, IL-4, IL-5, IL-6, IL-7, IL-9, I-TAC (CXCL11), KC (CXCL1), Leptin, Leptin R, LIX, L-Selectin (CD62L), Lungkine (CXCL15), Lymphotactin (XCL1), MCP-1 (CCL2), MCP-5, M-CSF, MDC (CCL22), MIG (CXCL9), MIP-1 alpha (CCL3), MIP-1 gamma, MIP-2, MIP-3 beta (CCL19), MIP-3 alpha (CCL20), MMP-2, MMP-3, Osteopontin (SPP1), Osteoprotegerin (TNFRSF11B), Platelet Factor 4 (CXCL4), Pro-MMP-9, P-Selectin, RANTES (CCL5), Resistin, SCF, SDF-1 alpha (CXCL12 alpha), Sonic Hedgehog N-Terminal (Shh-N), TNF RI (TNFRSF1A), TNF RII (TNFRSF1B), TARC (CCL17), I-309 (TCA-3/CCL1), TECK (CCL25), TCK-1 (CXCL7), TIMP-1, TIMP-2, TNF alpha, Thrombopoietin (TPO), TRANCE (TNFSF11), TROY (TNFRSF19), TSLP, VCAM-1 (CD106), VEGF-A, VEGFR1, VEGFR2, VEGFR3, VEGF-D

Cytokine arrays are an antibody-pair-based assay, analogous to ELISA, but using a membrane as a substrate rather than a plate. Capture antibodies are supplied arrayed/spotted on a membrane with each pair of spots representing a different analyte. Sample is added (0.2-1ml of 1 sample to each membrane), and then paired biotinylated detector antibodies and streptavidin HRP. The cytokine array is analyzed using the same methods as a chemiluminescent western blot. Comparison between samples can be by eye or using densitometry software for a semi-quantitative comparison.



## Images

---

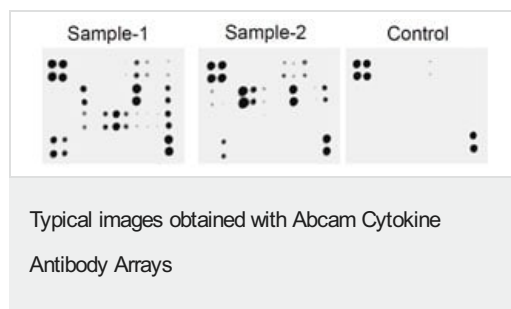

Typical images obtained with Abcam Mouse Cytokine Antibody Array - Membrane (96 Targets) ab193659. These membranes were probed with conditioned media from two different cell lines.

Membranes were exposed to film at room temperature for 1 minute.

**Please note:** All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

## Our Abpromise to you: Quality guaranteed and expert technical support

---

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
  
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <https://www.abcam.com/abpromise> or contact our technical team.

## Terms and conditions

---

- Guarantee only valid for products bought direct from Abcam or one of our authorized distributors