abcam

Product datasheet

Human NOS1 Antibody Pair - BSA and Azide free ab244100



3 Images

Overview

Product name Human NOS1 Antibody Pair - BSA and Azide free

Assay type ELISA set

312.5 pg/ml - 20000 pg/ml Range

Species reactivity Reacts with: Human

Product overview The Antibody Pair can be used to quantify Human NOS1. BSA and Azide free antibody pairs

> include unconjugated capture and detector antibodies suitable for sandwich ELISAs. The antibodies are provided at an approximate concentration of 1 mg/ml as measured by the protein A280 method. The recommended antibody orientation is based on internal optimization for ELISA-based assays. Antibody orientation is assay dependent and needs to be optimized for each assay type. Both capture and detector antibodies are rabbit monoclonal antibodies

delivering consistent, specific, and sensitive results.

For additional information on the performance of the antibody pair, see the equivalent SimpleStep ELISA® Kit (ab233628), which uses the same antibodies. However, due to differences in their formulation, this antibody pair cannot be used with the consumables provided with our SimpleStep ELISA Kits. Please note that the range provided for the pairs is only an estimation based on the performance of the related product using the same antibody pair. Performance of the antibody pair will depend on the specific characteristics of your assay. We guarantee the product works in sandwich ELISA, but we do not guarantee the sensitivity or dynamic range of the

antibody pair in your assay.

Tested applications Suitable for: Sandwich ELISA

Platform Reagents

Properties

Store at +4°C. Please refer to protocols. Storage instructions

Carrier free Yes

Components	10 x 96 tests
Human Niemann Pick C2 Capture Antibody (unconjugated)	1 x 100µg
Human Niemann Pick C2 Detector Antibody (unconjugated)	1 x 100µg

Function Produces nitric oxide (NO) which is a messenger molecule with diverse functions throughout the body. In the brain and peripheral nervous system, NO displays many properties of a neurotransmitter. Probably has nitrosylase activity and mediates cysteine S-nitrosylation of

cytoplasmic target proteins such SRR.

Tissue specificity Isoform 1 is ubiquitously expressed: detected in skeletal muscle and brain, also in testis, lung and

kidney, and at low levels in heart, adrenal gland and retina. Not detected in the platelets. Isoform 3 is expressed only in testis. Isoform 4 is detected in testis, skeletal muscle, lung, and kidney, at low

levels in the brain, but not in the heart and adrenal gland.

Sequence similarities Belongs to the NOS family.

Contains 1 FAD-binding FR-type domain.

Contains 1 flavodoxin-like domain.
Contains 1 PDZ (DHR) domain.

DomainThe PDZ domain in the N-terminal part of the neuronal isoform participates in protein-protein

interaction, and is responsible for targeting nNos to synaptic membranes in muscles. Mediates

interaction with VAC14.

Post-translational modifications

Ubiquitinated; mediated by STUB1/CHIP in the presence of Hsp70 and Hsp40 (in vitro).

Cellular localization Cell membrane > sarcolemma. Cell projection > dendritic spine. In skeletal muscle, it is localized

beneath the sarcolemma of fast-twitch muscle fiber by associating with the dystrophin

glycoprotein complex. In neurons, enriched in dendritic spines.


Applications

The Abpromise guarantee Our Abpromise guarantee covers the use of ab244100 in the following tested applications.

The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.


Application	Abreviews	Notes
Sandwich ELISA		Use at an assay dependent concentration.

Images

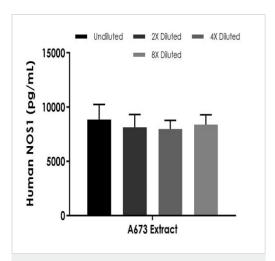


Sandwich ELISA - Human NOS1 Antibody Pair - BSA and Azide free (ab244100)

Representative standard curve from corresponding SimpleStep ELISA® Kit (ab233628), which uses the same antibody pair. For additional information on the performance of pair and kit, refer to the corresponding kit datasheet. Due to differences in the formulation and format of the antibodies in this pair, they cannot be used as substitutes for the antibody components in our SimpleStep ELISA® Kits.

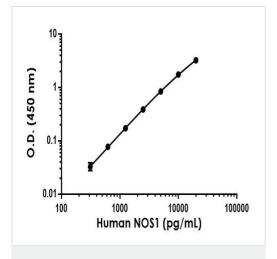

Sandwich ELISA - Human NOS1 Antibody Pair - BSA and Azide free (ab244100)

Representative standard curve from corresponding SimpleStep ELISA® Kit (ab233628), which uses the same antibody pair. For additional information on the performance of pair and kit, refer to the corresponding kit datasheet. Due to differences in the formulation and format of the antibodies in this pair, they cannot be used as substitutes for the antibody components in our SimpleStep ELISA® Kits.

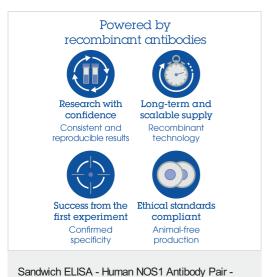

To learn more about the advantages of recombinant antibodies see here.

BSA and Azide free (ab244100)

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.com/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors