# abcam

# Product datasheet

# Mouse/Rat Retinol Binding Protein 4 Antibody Pair - BSA and Azide free ab241749



RabMAb

### 4 Images

#### Overview

**Product name** Mouse/Rat Retinol Binding Protein 4 Antibody Pair - BSA and Azide free

Assay type ELISA set

Range 15.62 pg/ml - 1000 pg/ml
Species reactivity Reacts with: Mouse, Rat

**Product overview** The Antibody Pair can be used to quantify Mouse/Rat RBP4. BSA and Azide free antibody pairs

include unconjugated capture and detector antibodies suitable for sandwich ELISAs. The antibodies are provided at an approximate concentration of 1 mg/ml as measured by the protein A280 method. The recommended antibody orientation is based on internal optimization for ELISA-based assays. Antibody orientation is assay dependent and needs to be optimized for each assay type. Both capture and detector antibodies are rabbit monoclonal antibodies

delivering consistent, specific, and sensitive results.

For additional information on the performance of the antibody pair, see the equivalent SimpleStep ELISA® Kit (ab202404), which uses the same antibodies. However, due to differences in their formulation, this antibody pair cannot be used with the consumables provided with our SimpleStep ELISA Kits. Please note that the range provided for the pairs is only an estimation based on the performance of the related product using the same antibody pair. Performance of the antibody pair will depend on the specific characteristics of your assay. We guarantee the product works in sandwich ELISA, but we do not guarantee the sensitivity or dynamic range of the

Tested applications Suitable for: Sandwich ELISA

**Platform** Reagents

**Properties** 

**Storage instructions** Store at +4°C. Please refer to protocols.

antibody pair in your assay.

Carrier free Yes

1

Components	10 x 96 tests
Mouse/Rat Retinol Binding Protein 4 Capture Antibody (unconjugated)	1 x 100µg
Mouse/Rat Retinol Binding Protein 4 Detector Antibody (unconjugated)	1 x 100µg

**Function**Delivers retinol from the liver stores to the peripheral tissues. In plasma, the RBP-retinol complex interacts with transthyretin, this prevents its loss by filtration through the kidney glomeruli.

**Involvement in disease** Defects in RBP4 are a cause of retinol-binding protein deficiency (RBP deficiency) [MIM:180250].

This condition causes night vision problems. It produces a typical 'fundus xerophthalmicus',

featuring a progressed atrophy of the retinal pigment epithelium.

**Sequence similarities**Belongs to the calycin superfamily. Lipocalin family.

Cellular localization Secreted.

#### **Applications**


The Abpromise guarantee

Our Abpromise guarantee covers the use of ab241749 in the following tested applications.


The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

Application	Abreviews	Notes
Sandwich ELISA		Use at an assay dependent concentration.

#### **Images**

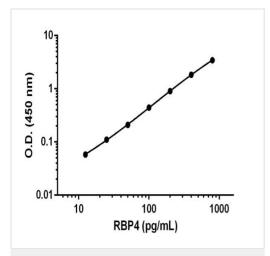


Sandwich ELISA - Mouse/Rat Retinol Binding Protein 4 Antibody Pair - BSA and Azide free (ab241749) Representative standard curve from corresponding SimpleStep ELISA® Kit (ab202404), which uses the same antibody pair. For additional information on the performance of pair and kit, refer to the corresponding kit datasheet. Due to differences in the formulation and format of the antibodies in this pair, they cannot be used as substitutes for the antibody components in our SimpleStep ELISA® Kits.

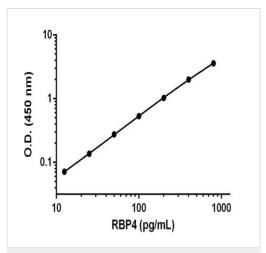


Sandwich ELISA - Mouse/Rat Retinol Binding Protein 4 Antibody Pair - BSA and Azide free (ab241749)

Representative standard curve from corresponding SimpleStep ELISA® Kit (ab202404), which uses the same antibody pair. For additional information on the performance of pair and kit, refer to the corresponding kit datasheet. Due to differences in the formulation and format of the antibodies in this pair, they cannot be used as substitutes for the antibody components in our SimpleStep ELISA® Kits.

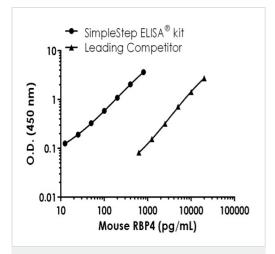

Sandwich ELISA - Mouse/Rat Retinol Binding Protein 4 Antibody Pair - BSA and Azide free (ab241749)

Representative standard curve from corresponding SimpleStep ELISA® Kit (ab202404), which uses the same antibody pair. For additional information on the performance of pair and kit, refer to the corresponding kit datasheet. Due to differences in the formulation and format of the antibodies in this pair, they cannot be used as substitutes for the antibody components in our SimpleStep ELISA® Kits.

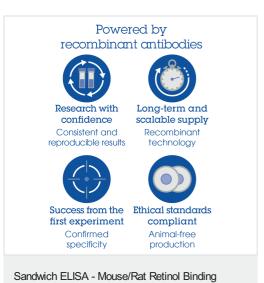

Protein 4 Antibody Pair - BSA and Azide free

(ab241749)

To learn more about the advantages of recombinant antibodies see **here**.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

#### Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit <a href="https://www.abcam.com/abpromise">https://www.abcam.com/abpromise</a> or contact our technical team.

## Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors