abcam

Product datasheet

Alexa Fluor® 647 Anti-Hemopexin antibody [EPR5610] ab197882



3 Images

Overview

Product name Alexa Fluor® 647 Anti-Hemopexin antibody [EPR5610]

Description Alexa Fluor® 647 Rabbit monoclonal [EPR5610] to Hemopexin

Host species Rabbit

Conjugation Alexa Fluor® 647. Ex: 652nm. Em: 668nm

Tested applications Suitable for: ICC/IF, Flow Cyt (Intra)

Species reactivity Reacts with: Human

Does not react with: Mouse. Rat

Synthetic peptide. This information is proprietary to Abcam and/or its suppliers. **Immunogen**

Positive control Flow Cyt (intra): HepG2 cells. ICC/IF: HepG2 cells.

General notes This product is a recombinant monoclonal antibody, which offers several advantages including:

- High batch-to-batch consistency and reproducibility

- Improved sensitivity and specificity

- Long-term security of supply

- Animal-free production

For more information see here.

Our RabMAb® technology is a patented hybridoma-based technology for making rabbit monoclonal antibodies. For details on our patents, please refer to **RabMAb**® **patents**.


Properties

Form Liquid

Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C. Storage instructions

Avoid freeze / thaw cycle. Store In the Dark.

 $K_D = 8.30 \times 10^{-12} M$ Dissociation constant (K_D)

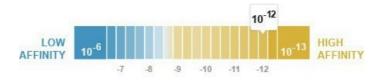

Learn more about K_D

Storage buffer pH: 7.40

Preservative: 0.02% Sodium azide

Constituents: 30% Glycerol (glycerin, glycerine), 1% BSA, PBS

Purity Protein A purified

Clonality Monoclonal
Clone number EPR5610

Isotype IgG

Applications


The Abpromise guarantee Our Abpromise guarantee covers the use of ab197882 in the following tested applications.

The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

Application	Abreviews	Notes
ICC/IF		1/100. This product gave a positive signal in HepG2 cells fixed with 4% formaldehyde (10 min) and 100% methanol (5 min).
Flow Cyt (Intra)		1/50. ab199093 - Rabbit monoclonal lgG (Alexa Fluor® 647), is suitable for use as an isotype control with this antibody.

Target		
Function	Binds heme and transports it to the liver for breakdown and iron recovery, after which the free hemopexin returns to the circulation.	
Tissue specificity	Expressed by the liver and secreted in plasma.	
Sequence similarities	Belongs to the hemopexin family. Contains 5 hemopexin-like domains.	
Post-translational modifications	N- and O-glycosylated. O-glycosylated with core 1 or possibly core 8 glycans.	
Cellular localization	Secreted.	

Images

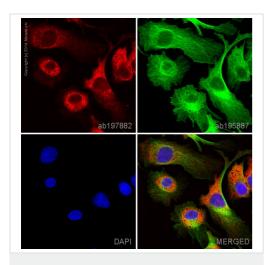


Immunocytochemistry/ Immunofluorescence - Alexa Fluor® 647 Anti-Hemopexin antibody [EPR5610] (ab197882)

ab197882 staining Hemopexin in HepG2 cells. The cells were fixed with 4% formaldehyde (10 min), permeabilized with 0.1% Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1% PBS-Tween for 1h. The cells were then incubated overnight at +4°C with ab197882 at a 1/100 dilution (shown in red) and ab195887, Mouse monoclonal to alpha Tubulin (Alexa Fluor[®] 488), at a 1/250 dilution (shown in green). Nuclear DNA was labelled with DAPI (shown in blue).

Image was taken with a confocal microscope (Leica-Microsystems, TCS SP8).

This product also gave a positive signal under the same testing conditions in HepG2 cells fixed with 100% methanol (5 min)

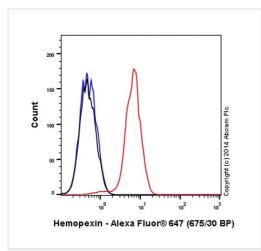


Flow Cytometry (Intracellular) - Alexa Fluor® 647 Anti-Hemopexin antibody [EPR5610] (ab197882)

Overlay histogram showing HepG2 cells stained with ab197882 (red line). The cells were fixed with 4% formaldehyde (10 min) and then permeabilized with 0.1% PBS-Tween for 20 min. The cells were then incubated in 1x PBS / 10% normal goat serum / 0.3M glycine to block non-specific protein-protein interactions followed by the antibody (ab197882, 1/50 dilution) for 30 min at 22°C. Isotype control antibody (black line) was rabbit IgG (monoclonal) Alexa Fluor[®] 647 used at the same concentration and conditions as the primary antibody. Unlabelled sample (blue line) was also used as a control.

Acquisition of >5,000 events were collected using a solid-state 25mW red diode laser (635 nm) and 675/30 bandpass filter.

This antibody gave a positive signal in HepG2 fixed with 80% methanol (5 min)/permeabilized with 0.1% PBS-Tween for 20 min used under the same conditions.

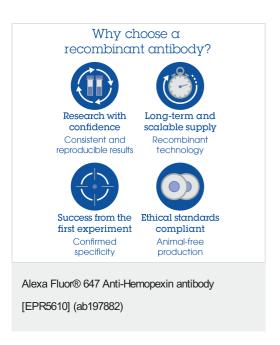

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

Our Abpromise to you: Quality guaranteed and expert technical support

- Replacement or refund for products not performing as stated on the datasheet
- Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.com/abpromise or contact our technical team.

Terms and conditions

• Guarantee only valid for products bought direct from Abcam or one of our authorized distributors