abcam

Product datasheet

Anti-NuMA antibody ab86129

*** * * 2 Abreviews 2 References 3 Images

Overview

Product name Anti-NuMA antibody

Description Rabbit polyclonal to NuMA

Host species Rabbit

Tested applications Suitable for: ICC, WB, IHC-P

Species reactivity Reacts with: Human

Immunogen Synthetic peptide. This information is proprietary to Abcam and/or its suppliers.

Positive control WB: MCF7 whole cell lysate. ICC: HeLa cells. IHC: human breast adenocarcinoma.

General notes

The Life Science industry has been in the grips of a reproducibility crisis for a number of years.

Abcam is leading the way in addressing this with our range of recombinant monoclonal antibodies and knockout edited cell lines for gold-standard validation. Please check that this product meets

your needs before purchasing.

If you have any questions, special requirements or concerns, please send us an inquiry and/or contact our Support team ahead of purchase. Recommended alternatives for this product can be

found below, along with publications, customer reviews and Q&As

Properties

Form Liquid

Storage instructions Shipped at 4°C. Store at +4°C short term (1-2 weeks). Upon delivery aliquot. Store at -20°C or -

80°C. Avoid freeze / thaw cycle.

Storage buffer pH: 7.40

Preservative: 0.02% Sodium azide

Constituent: PBS

Batches of this product that have a concentration < 1mg/ml may have BSA added as a stabilising

agent. If you would like information about the formulation of a specific lot, please contact our

scientific support team who will be happy to help.

Purity Immunogen affinity purified

Clonality Polyclonal

Isotype IgG

Applications

The Abpromise guarantee

Our Abpromise guarantee covers the use of ab86129 in the following tested applications.

The application notes include recommended starting dilutions; optimal dilutions/concentrations should be determined by the end user.

Application	Abreviews	Notes
ICC		Use a concentration of 1 µg/ml.
WB		Use a concentration of 1 µg/ml. Detects a band of approximately 235 kDa (predicted molecular weight: 238 kDa).
IHC-P	**** <u>(1)</u>	Use a concentration of 1 µg/ml. Perform heat mediated antigen retrieval before commencing with IHC staining protocol.

Target


Function

Cellular localization

May be a structural component of the nucleus.

Nucleus. Chromosome. Dissociates from condensing chromosomes during early prophase, before the complete disintegration of the nuclear lamina. As mitosis progresses it reassociates with telophase chromosomes very early during nuclear reformation, before substantial accumulation of lamins on chromosomal surfaces is evident.


Images

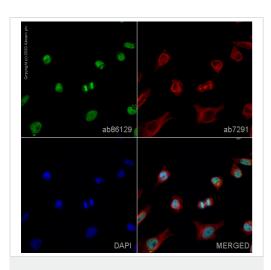

Immunocytochemistry - Anti-NuMA antibody (ab86129)

ab86129 staining NuMa in HeLa cells. The cells were fixed with 4% paraformaldehyde (10 min), permeabilized with 0.1% PBS-Triton X-100 for 5 minutes and then blocked with 1% BSA/10% normal goat serum/0.3M glycine in 0.1%PBS-Tween for 1h. The cells were then incubated overnight at 4°C with ab86129 at 1µg/ml and ab7291, Mouse monoclonal [DM1A] to alpha Tubulin - Loading Control. Cells were then incubated with ab150081, Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (Alexa Fluor® 488), pre-adsorbed at 1/1000 dilution (shown in green) and ab150080, Goat polyclonal Secondary Antibody to Rabbit IgG - H&L (Alexa Fluor® 594) at 1/1000 dilution (shown in pseudocolour red). Nuclear DNA was labelled with DAPI (shown in blue).

Image was acquired with a high-content analyser (Operetta CLS, Perkin Elmer) and a maximum intensity projection of confocal sections is shown.

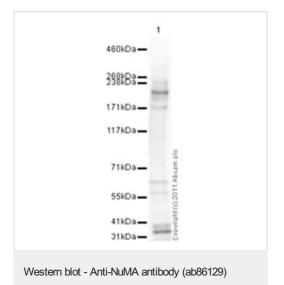

Anti-NuMA antibody (ab86129) at 1 µg/ml + MCF7 (Human breast adenocarcinoma cell line) Whole Cell Lysate at 10 µg

Secondary

Goat Anti-Rabbit IgG H&L (HRP) preadsorbed (ab97080) at 1/5000 dilution

Developed using the ECL technique.

Performed under reducing conditions.

Predicted band size: 238 kDa Observed band size: 235 kDa


Additional bands at: 171 kDa, 35 kDa, 65 kDa. We are unsure as

to the identity of these extra bands.

Exposure time: 16 minutes

Abcam recommends using a transfer buffer including 20% Ethanol and SDS as well as unheating tissue lysates for this product.

Abcam welcomes customer feedback and would appreciate any comments regarding this product and the data presented above.

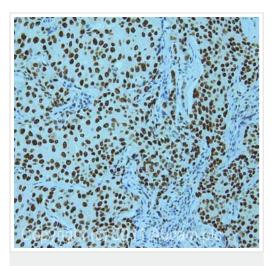

Immunohistochemistry (Formalin/PFA-fixed paraffinembedded sections) - Anti-NuMA antibody (ab86129)

IHC image of NuMA staining in human breast adenocarcinoma formalin fixed paraffin embedded tissue section, performed on a Leica Bond TM system using the standard protocol F. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH6, epitope retrieval solution 1) for 20 mins. The section was then incubated with ab86129, 1µg/ml, for 15 mins at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX. For other IHC staining systems (automated and non-automated) customers should optimize variable parameters such as antigen retrieval conditions, primary antibody concentration and antibody incubation times

Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- · Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- · We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.com/abpromise or contact our technical team.

Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors