

Version 1 Last updated 11 September 2017

# ab223862 Mouse HGF SimpleStep ELISA<sup>®</sup> Kit

For the quantitative measurement of HGF in Mouse serum, plasma, cell culture supernatant, cell and tissue extract samples.

This product is for research use only and is not intended for diagnostic use.

## Table of Contents

1. Overview	1
2. Protocol Summary	2
3. Precautions	3
4. Storage and Stability	3
5. Limitations	4
6. Materials Supplied	4
7. Materials Required, Not Supplied	5
8. Technical Hints	5
9. Reagent Preparation	7
10. Standard Preparation	8
11. Sample Preparation	9
12. Plate Preparation	12
13. Assay Procedure	13
14. Calculations	15
15. Typical Data	16
16. Typical Sample Values	18
17. Assay Specificity	28
18. Species Reactivity	29
19. Troubleshooting	30
20. Notes	31
<b>Technical Support</b>	<b>34</b>

## 1. Overview

HGF *in vitro* SimpleStep ELISA<sup>®</sup> (Enzyme-Linked Immunosorbent Assay) kit is designed for the quantitative measurement of HGF protein in Mouse serum, plasma, cell culture supernatant, cell and tissue extract samples.

The SimpleStep ELISA<sup>®</sup> employs an affinity tag labeled capture antibody and a reporter conjugated detector antibody which immunocapture the sample analyte in solution. This entire complex (capture antibody/analyte/detector antibody) is in turn immobilized via immunoaffinity of an anti-tag antibody coating the well. To perform the assay, samples or standards are added to the wells, followed by the antibody mix. After incubation, the wells are washed to remove unbound material. TMB substrate is added and during incubation is catalyzed by HRP, generating blue coloration. This reaction is then stopped by addition of Stop Solution completing any color change from blue to yellow. Signal is generated proportionally to the amount of bound analyte and the intensity is measured at 450 nm. Optionally, instead of the endpoint reading, development of TMB can be recorded kinetically at 600 nm.

Hepatocyte Growth Factor (HGF) is mitogen of parenchymal hepatocytes. HGF also acts as a growth factor and is expressed by a wide variety of tissue and cell types. HGF binds to the HGF-Receptor to help regulate cell growth, cell motility and morphogenesis. Elevated circulating levels of HGF can be measured as prognostic marker of severity of hypertension and heart failure. HGF is secreted as a full length inactive polypeptide, and then is cleaved by serine proteases into a 69-kDa alpha-chain and 34-kDa beta chain which dimerize producing the active heterodimeric molecule. Mouse HGF SimpleStep ELISA<sup>®</sup> (ab223862) is designed to recognize this dimerized protein.

## 2. Protocol Summary

Prepare all reagents, samples, and standards as instructed



Add 50  $\mu$ L standard or sample to appropriate wells



Add 50  $\mu$ L Antibody Cocktail to all wells



Incubate at room temperature for 1 hour



Aspirate and wash each well three times with 350  $\mu$ L 1X Wash Buffer  
PT



Add 100  $\mu$ L TMB Substrate to each well and incubate for 10 minutes.



Add 100  $\mu$ L Stop Solution and read OD at 450 nm

### 3. Precautions

**Please read these instructions carefully prior to beginning the assay.**

- All kit components have been formulated and quality control tested to function successfully as a kit.
- We understand that, occasionally, experimental protocols might need to be modified to meet unique experimental circumstances. However, we cannot guarantee the performance of the product outside the conditions detailed in this protocol booklet.
- Reagents should be treated as possible mutagens and should be handled with care and disposed of properly. Please review the Safety Datasheet (SDS) provided with the product for information on the specific components.
- Observe good laboratory practices. Gloves, lab coat, and protective eyewear should always be worn. Never pipet by mouth. Do not eat, drink or smoke in the laboratory areas.
- All biological materials should be treated as potentially hazardous and handled as such. They should be disposed of in accordance with established safety procedures.

### 4. Storage and Stability

**Store kit at +4°C immediately upon receipt. Kit has a storage time of 1 year from receipt, providing components have not been reconstituted.**

Refer to list of materials supplied for storage conditions of individual components.

## 5. Limitations

- Assay kit intended for research use only. Not for use in diagnostic procedures.
- Do not mix or substitute reagents or materials from other kit lots or vendors. Kits are QC tested as a set of components and performance cannot be guaranteed if utilized separately or substituted.

## 6. Materials Supplied

Item	Quantity	Storage Condition
Mouse HGF Capture Antibody 10X	600 µL	+4°C
Mouse HGF Detector Antibody 10X	600 µL	+4°C
Mouse HGF Lyophilized Recombinant Protein	2 Vials	+4°C
Antibody Diluent CPR	6 mL	+4°C
Wash Buffer PT 10X	20 mL	+4°C
Cell Extraction Buffer PTR 5X	10 mL	+4°C
Cell Extraction Enhancer Solution 50X	1 mL	+4°C
TMB Substrate	12 mL	+4°C
Stop Solution	12 mL	+4°C
Sample Diluent NS	50 mL	+4°C
Anti-tag coated microplate (12 x 8 well strips)	96 Wells	+4°C
Plate Seal	1	+4°C

## 7. Materials Required, Not Supplied

These materials are not included in the kit, but will be required to successfully perform this assay:

- Microplate reader capable of measuring absorbance at 450 or 600 nm.
- Method for determining protein concentration (BCA assay recommended).
- Deionized water.
- Multi- and single-channel pipettes.
- Tubes for standard dilution.
- Plate shaker for all incubation steps.
- Optional: Phenylmethylsulfonyl Fluoride (PMSF) (or other protease inhibitors).

## 8. Technical Hints

- Samples generating values higher than the highest standard should be further diluted in the appropriate sample dilution buffers.
- Avoid foaming or bubbles when mixing or reconstituting components.
- Avoid cross contamination of samples or reagents by changing tips between sample, standard and reagent additions.
- Ensure plates are properly sealed or covered during incubation steps.
- Complete removal of all solutions and buffers during wash steps is necessary to minimize background.
- As a guide, typical ranges of sample concentration for commonly used sample types are shown below in Sample Preparation (section 11).
- All samples should be mixed thoroughly and gently.
- Avoid multiple freeze/thaw of samples.
- Incubate ELISA plates on a plate shaker during all incubation steps.
- When generating positive control samples, it is advisable to change pipette tips after each step.

- The provided Antibody Diluents and Sample Diluents contain protease inhibitor aprotinin. Additional protease inhibitors can be added if required.
- The provided Cell Extraction Buffer 5X contains phosphatase inhibitors and protease inhibitor aprotinin. Additional protease inhibitors can be added if required.
- The provided Cell Extraction Enhancer Solution 50X may precipitate when stored at + 4°C. To dissolve, warm briefly at + 37°C and mix gently. The Cell Extraction Enhancer Solution 50X can be stored at room temperature to avoid precipitation.
- **To avoid high background always add samples or standards to the well before the addition of the antibody cocktail.**
- **This kit is sold based on number of tests. A 'test' simply refers to a single assay well. The number of wells that contain sample, control or standard will vary by product. Review the protocol completely to confirm this kit meets your requirements. Please contact our Technical Support staff with any questions.**

## 9. Reagent Preparation

- Equilibrate all reagents to room temperature (18-25°C) prior to use. The kit contains enough reagents for 96 wells. The sample volumes below are sufficient for 48 wells (6 x 8-well strips); adjust volumes as needed for the number of strips in your experiment.
- Prepare only as much reagent as is needed on the day of the experiment. Capture and Detector Antibodies have only been tested for stability in the provided 10X formulations.

### 9.1 1X Cell Extraction Buffer PTR (For cell and tissue extracts only):

Prepare 1X Cell Extraction Buffer PTR by diluting Cell Extraction Buffer PTR 5X and 50X Cell Extraction Enhancer Solution to 1X with deionized water. To make 10 mL 1X Cell Extraction Buffer PTR combine 7.8 mL deionized water, 2 mL Cell Extraction Buffer PTR 5X and 200  $\mu$ L Cell Extraction Enhancer Solution 50X. Mix thoroughly and gently. If required protease inhibitors can be added.

Alternative – Enhancer may be added to 1X Cell Extraction Buffer PTR after extraction of cells or tissue. Refer to note in the Troubleshooting section.

### 9.2 1X Wash Buffer PT:

Prepare 1X Wash Buffer PT by diluting Wash Buffer PT 10X with deionized water. To make 50 mL 1X Wash Buffer PT combine 5 mL Wash Buffer PT 10X with 45 mL deionized water. Mix thoroughly and gently.

### 9.3 Antibody Cocktail:

Prepare Antibody Cocktail by diluting the capture and detector antibodies in Antibody Diluent CPR. To make 3 mL of the Antibody Cocktail combine 300  $\mu$ L 10X Capture Antibody and 300  $\mu$ L 10X Detector Antibody with 2.4 mL Antibody Diluent CPR. Mix thoroughly and gently.

## 10. Standard Preparation

- Always prepare a fresh set of standards for every use.
- Discard working standard dilutions after use as they do not store well.
- The following section describes the preparation of a standard curve for duplicate measurements (recommended).

**10.1 For serum, plasma and cell culture supernatant sample measurements,** reconstitute the HGF protein standard by adding 500  $\mu\text{L}$  of Sample Diluent NS.


For **cell and tissue sample measurements,** reconstitute the HGF protein standard by adding 500  $\mu\text{L}$  of 1X Cell Extraction Buffer PTR.

Hold at room temperature for 10 minutes and mix thoroughly and gently. This is the 10,000  $\text{pg/mL}$  **Stock Standard Solution**.

**10.2** Label eight tubes, Standards 1– 8.

**10.3** Add 180  $\mu\text{L}$  of appropriate diluent (see step 10.1) into tube number 1 and 150  $\mu\text{L}$  of appropriate diluent into numbers 2-8.

**10.4** Use the Stock Standard to prepare the following dilution series. Standard #8 contains no protein and is the Blank control:

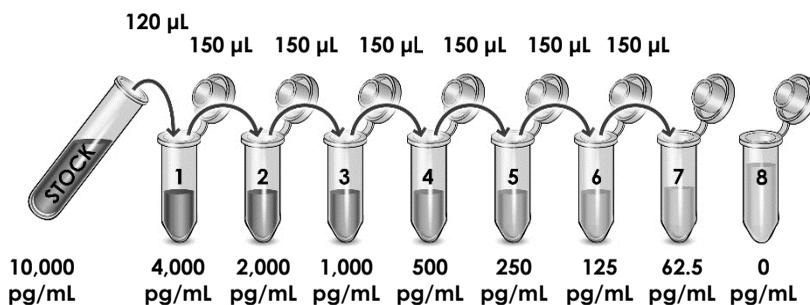

## 11. Sample Preparation

Typical Sample Dynamic Range	
Sample Type	Range
Mouse Plasma - Heparin	6 - 50%
Mouse Plasma - EDTA	3 - 50%
Mouse Plasma - Citrate	3 - 50%
Mouse Serum	1 - 20%
Cell Culture Media	50%
NIH/3T3 Cell Extract	125 – 2,000 µg/mL
3T3-L1 Cell Extract	62.5 – 1,000 µg/mL
Mouse Kidney Tissue Extract	1,000 – 4,000 µg/mL
Mouse Spleen Tissue Extract	250 – 4,000 µg/mL
Rat Plasma-Heparin	6 – 100%
Rat Plasma-EDTA	6 – 100%
Rat Plasma – Citrate	6 – 100%
Rat Serum	1 - 20%

### 11.1 Plasma:

Collect plasma using citrate, EDTA or heparin. Centrifuge samples at 2,000 x g for 10 minutes. Dilute samples into Sample Diluent NS and assay. Store un-diluted plasma samples at -20°C or below for up to 3 months. Avoid repeated freeze-thaw cycles.

### 11.2 Serum:

Samples should be collected into a serum separator tube. After clot formation, centrifuge samples at 2,000 x g for 10 minutes and collect serum. Dilute samples into Sample Diluent NS and assay. Store un-diluted serum at -20°C or below. Avoid repeated freeze-thaw cycles.

### 11.3 Cell Culture Supernatants:

Centrifuge cell culture media at 2,000 x g for 10 minutes to remove debris. Collect supernatants and assay. Or dilute samples into Sample Diluent NS and assay. Store un-diluted samples at -20°C or below. Avoid repeated freeze-thaw cycles.

### 11.4 Preparation of extracts from cell pellets:

- 11.4.1 Collect non-adherent cells by centrifugation or scrape to collect adherent cells from the culture flask. Typical centrifugation conditions for cells are 500 x g for 5 minutes at 4°C.
- 11.4.2 Rinse cells twice with PBS.
- 11.4.3 Solubilize pellet at  $2 \times 10^7$  cell/mL in chilled 1X Cell Extraction Buffer PTR.
- 11.4.4 Incubate on ice for 20 minutes.
- 11.4.5 Centrifuge at 18,000 x g for 20 minutes at 4°C.
- 11.4.6 Transfer the supernatants into clean tubes and discard the pellets.
- 11.4.7 Assay samples immediately or aliquot and store at -80°C. The sample protein concentration in the extract may be quantified using a protein assay.
- 11.4.8 Dilute samples to desired concentration in 1X Cell Extraction Buffer PTR.

### 11.5 Preparation of extracts from adherent cells by direct lysis (alternative protocol):

- 11.5.1 Remove growth media and rinse adherent cells 2 times in PBS.
- 11.5.2 Solubilize the cells by addition of chilled 1X Cell Extraction Buffer PTR directly to the plate (use 750  $\mu$ L - 1.5 mL 1X Cell Extraction Buffer PTR per confluent 15 cm diameter plate).
- 11.5.3 Scrape the cells into a microfuge tube and incubate the lysate on ice for 15 minutes.
- 11.5.4 Centrifuge at 18,000 x g for 20 minutes at 4°C.
- 11.5.5 Transfer the supernatants into clean tubes and discard the pellets.
- 11.5.6 Assay samples immediately or aliquot and store at -80°C. The sample protein concentration in the extract may be quantified using a protein assay.
- 11.5.7 Dilute samples to desired concentration in 1X Cell Extraction Buffer PTR.

## **11.6 Preparation of extracts from tissue homogenates:**

- 11.6.1 Tissue lysates are typically prepared by homogenization of tissue that is first minced and thoroughly rinsed in PBS to remove blood (dounce homogenizer recommended).
- 11.6.2 Homogenize 100 to 200 mg of wet tissue in 500  $\mu$ L – 1 mL of chilled 1X Cell Extraction Buffer PTR. For lower amounts of tissue adjust volumes accordingly.
- 11.6.3 Incubate on ice for 20 minutes.
- 11.6.4 Centrifuge at 18,000 x g for 20 minutes at 4°C.
- 11.6.5 Transfer the supernatants into clean tubes and discard the pellets.
- 11.6.6 Assay samples immediately or aliquot and store at -80°C. The sample protein concentration in the extract may be quantified using a protein assay.
- 11.6.7 Dilute samples to desired concentration in 1X Cell Extraction Buffer PTR.

## 12. Plate Preparation

- The 96 well plate strips included with this kit are supplied ready to use. It is not necessary to rinse the plate prior to adding reagents.
- Unused plate strips should be immediately returned to the foil pouch containing the desiccant pack, resealed and stored at 4°C.
- For each assay performed, a minimum of two wells must be used as the zero control.
- For statistical reasons, we recommend each sample should be assayed with a minimum of two replicates (duplicates).
- Differences in well absorbance or “edge effects” have not been observed with this assay.

## 13. Assay Procedure

- Equilibrate all materials and prepared reagents to room temperature prior to use.
- We recommend that you assay all standards, controls and samples in duplicate.
- Prepare all reagents, working standards, and samples as directed in the previous sections.

**13.1** Prepare all reagents, working standards, and samples as directed in the previous sections.

**13.2** Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, reseal and return to 4°C storage.

**13.3** Add 50 µL of all sample or standard to appropriate wells.

**13.4** Add 50 µL of the Antibody Cocktail to each well.

**13.5** Seal the plate and incubate for 1 hour at room temperature on a plate shaker set to 400 rpm.

**13.6** Wash each well with 3 x 350 µL 1X Wash Buffer PT. Wash by aspirating or decanting from wells then dispensing 350 µL 1X Wash Buffer PT into each well. Complete removal of liquid at each step is essential for good performance. After the last wash invert the plate and blot it against clean paper towels to remove excess liquid.

**13.7** Add 100 µL of TMB Substrate to each well and incubate for 10 minutes in the dark on a plate shaker set to 400 rpm.

**13.8** Add 100 µL of Stop Solution to each well. Shake plate on a plate shaker for 1 minute to mix. Record the OD at 450 nm. This is an endpoint reading.

**13.9** Alternative to 13.7 – 13.8: Instead of the endpoint reading at 450 nm, record the development of TMB Substrate kinetically. Immediately after addition of TMB Development Solution begin recording the blue color development with elapsed time in the microplate reader prepared with the following settings:

Mode	Kinetic
Wavelength:	600 nm
Time:	up to 15 min
Interval:	20 sec - 1 min
Shaking:	Shake between readings

**Δ Note:** that an endpoint reading can also be recorded at the completion of the kinetic read by adding 100  $\mu$ L Stop Solution to each well and recording the OD at 450 nm.


**13.10** Analyze the data as described below.

## 14. Calculations

- 14.1 Calculate the average absorbance value for the blank control (zero) standards. Subtract the average blank control standard absorbance value from all other absorbance values.
- 14.2 **Create a standard curve** by plotting the average blank control subtracted absorbance value for each standard concentration (y-axis) against the target protein concentration (x-axis) of the standard. Use graphing software to draw the best smooth curve through these points to construct the standard curve.
  - Δ **Note:** Most microplate reader software or graphing software will plot these values and fit a curve to the data. A four parameter curve fit (4PL) is often the best choice; however, other algorithms (e.g. linear, semi-log, log/log, 4 parameter logistic) can also be tested to determine if it provides a better curve fit to the standard values.
- 14.3 Determine the concentration of the target protein in the sample by interpolating the blank control subtracted **absorbance values against the standard curve**. Multiply the resulting value by the appropriate sample dilution factor, if used, to obtain the concentration of target protein in the sample.
- 14.4 Samples generating absorbance values greater than that of the highest standard should be further diluted and reanalyzed. Similarly, samples which measure at an absorbance values less than that of the lowest standard should be retested in a less dilute form.


## 15. Typical Data

Typical standard curve – data provided for demonstration purposes only. A new standard curve must be generated for each assay performed.

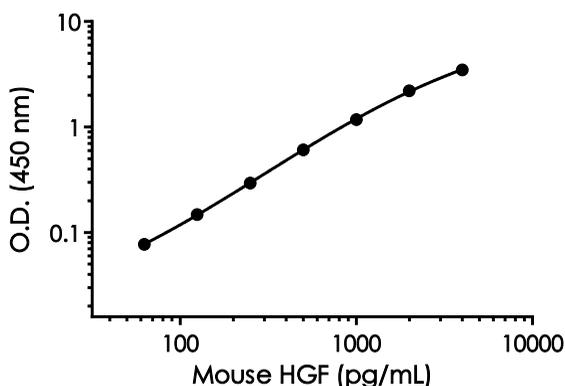

Standard Curve Measurements			
Concentration (pg/mL)	O.D 450 nm		Mean O.D
	1	2	
0	0.051	0.051	0.051
62.5	0.127	0.130	0.128
125	0.198	0.201	0.199
250	0.345	0.346	0.345
500	0.651	0.669	0.660
1000	1.216	1.241	1.229
2000	2.292	2.220	2.256
4000	3.507	3.568	3.537

**Figure 1.** Example of MouseHGF standard curve in Sample Diluent NS. The HGF standard curve was prepared as described in Section 10. Raw data values are shown in the table. Background-subtracted data values (mean +/- SD) are graphed.

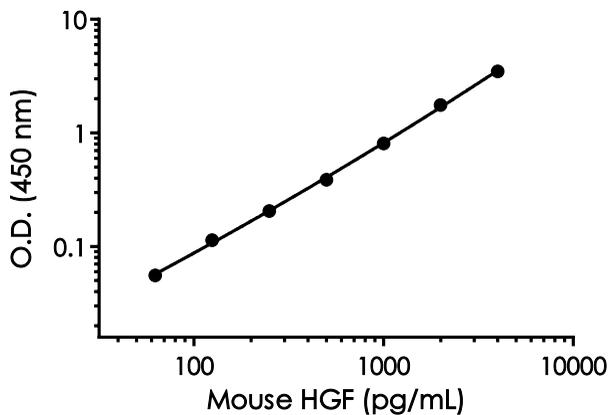

Standard Curve Measurements			
Concentration (pg/mL)	O.D 450 nm		Mean O.D
	1	2	
0	0.044	0.043	0.043
62.5	0.097	0.100	0.099
125	0.164	0.150	0.157
250	0.253	0.244	0.249
500	0.428	0.434	0.431
1000	0.856	0.846	0.851
2000	1.820	1.787	1.803
4000	3.511	3.534	3.522

**Figure 2.** Example of Mouse HGF standard curve in 1X Cell Extraction Buffer PTR. The HGF standard curve was prepared as described in Section 10. Raw data values are shown in the table. Background-subtracted data values (mean +/- SD) are graphed.

## 16. Typical Sample Values

### **SENSITIVITY –**

The calculated minimal detectable dose (MDD) is 8.4 pg/mL in Sample Diluent NS. The MDD was determined by calculating the mean of zero standard replicates (n=33) and adding 2 standard deviations then extrapolating the corresponding concentration.

The calculated minimal detectable dose (MDD) is 11 pg/mL in 1X Cell Extraction Buffer PTR. The MDD was determined by calculating the mean of zero standard replicates (n=6) and adding 2 standard deviations then extrapolating the corresponding concentration.

## RECOVERY –

Three concentrations of HGF Recombinant Protein were spiked in duplicate to the indicated biological matrix to evaluate signal recovery in the working range of the assay.

<b>Sample Type</b>	<b>Average % Recovery</b>	<b>Range (%)</b>
5% Mouse Serum	111	109-117
50% Mouse Plasma-EDTA	104	100-108
50% Mouse Plasma – Citrate	108	105-111
50% Mouse Plasma-Heparin	96	93-98
500 µg/mL Mouse Kidney Tissue Extract	119	113-123
500 µg/mL Mouse Spleen Tissue Extract	119	115-122
250 µg/mL NIH/3T3 Cell Extract	104	101-107
50% Cell Culture Media	109	107-112

## Linearity of Dilution

Linearity of dilution is determined based on interpolated values from the standard curve. Linearity of dilution defines a sample concentration interval in which interpolated target concentrations are directly proportional to sample dilution.

Native HGF was measured in the following biological samples in a 2-fold dilution series. Sample dilutions are made in Sample Diluent NS.

Dilution Factor	Interpolated value	20% Mouse Serum	50% Mouse Plasma (Citrate)	50% Mouse Plasma (EDTA)	50% Mouse Plasma (Heparin)
Undiluted	pg/mL	2524	306.0	527.2	416.9
	<b>% Expected value</b>	100	100	100	100
2	pg/mL	1114	185.1	286.7	241.8
	<b>% Expected value</b>	88	121	109	116
4	pg/mL	586	88.5	154.9	108.9
	<b>% Expected value</b>	93	116	118	104
8	pg/mL	269	44.3	69.3	56.9
	<b>% Expected value</b>	91	116	105	109
16	pg/mL	136	21.4	29.6	NL
	<b>% Expected value</b>	86	112	90	-

NL – Non linear

Native HGF was measured in the following biological samples in a 2-fold dilution series. Sample dilutions are made in Sample Diluent NS.

Dilution Factor	Interpolated value	20% Rat Serum	100% Rat Plasma (Citrate)	100% Rat Plasma (EDTA)	100% Rat Plasma (Heparin)
Undiluted	pg/mL	3829	2401	3406	2729
	<b>% Expected value</b>	100	100	100	100
2	pg/mL	2297	1171	1607	1289
	<b>% Expected value</b>	120	98	94	94
4	pg/mL	946	601	466	666
	<b>% Expected value</b>	99	100	90	98
8	pg/mL	466	313	385	322
	<b>% Expected value</b>	97	104	90	94
16	pg/mL	247	162	208	169
	<b>% Expected value</b>	103	108	98	99

Native HGF was measured in the following biological samples in a 2-fold dilution series. Sample dilutions are made in 1X Cell Extraction Buffer PTR.

Dilution Factor	Interpolated value	2000 µg/mL NIH/3T3 Cell Extract	1000 µg/mL 3T3-L1 Cell Extract	4000 µg/mL Mouse Spleen Tissue Extract	4000 µg/mL Mouse Kidney Tissue Extract
Undiluted	pg/mL	4053	2972	2568	1386
	<b>% Expected value</b>	100	100	100	100
2	pg/mL	2184	1349	1353	659
	<b>% Expected value</b>	108	91	105	95
4	pg/mL	1022	638	646	293
	<b>% Expected value</b>	101	86	101	84
8	pg/mL	459	309	283	NL
	<b>% Expected value</b>	91	83	88	-
16	pg/mL	217	156	140	-
	<b>% Expected value</b>	86	84	87	-

NL – Non linear


Recombinant HGF was spiked into the following biological samples and diluted in a 2-fold dilution series in Sample Diluent NS.

Dilution Factor	Interpolated value	50% Cell Culture Media
Undiluted	pg/mL	1984
	<b>% Expected value</b>	100
2	pg/mL	933
	<b>% Expected value</b>	94
4	pg/mL	466
	<b>% Expected value</b>	94
8	pg/mL	242
	<b>% Expected value</b>	98
16	pg/mL	120
	<b>% Expected value</b>	97

#### PRECISION –


Mean coefficient of variations of interpolated values from mouse serum within the working range of the assay.

	Intra-Assay	Inter-Assay
n =	5	3
CV(%)	3.2	4.9

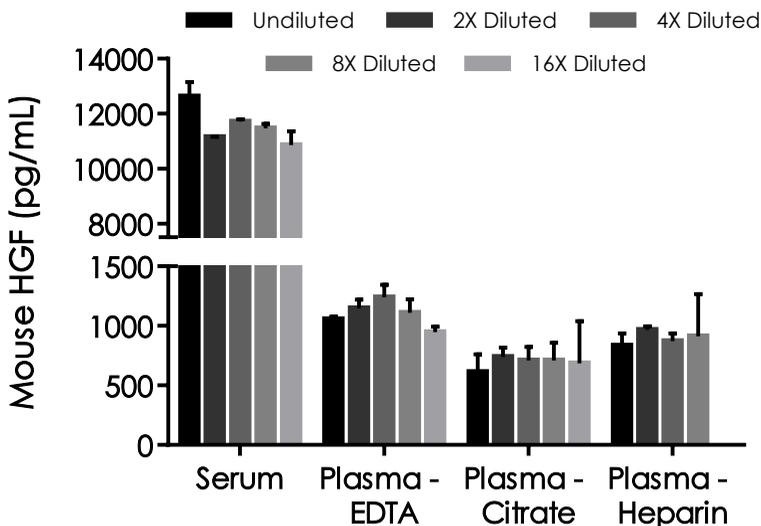


Fi


**Figure 3.** Interpolated concentrations of native HGF in Mouse serum and plasma samples. The concentrations of HGF were measured in duplicates, interpolated from the HGF standard curves and corrected for sample dilution. Undiluted samples are as follows: serum 20%, plasma (EDTA) 50%, plasma (citrate) 50%, and plasma (heparin) 50%. The interpolated dilution factor corrected values are plotted (mean +/- SD, n=2). The mean HGF concentration was determined to be 11,764 pg/mL in serum, 1099 pg/mL in plasma (EDTA), 691 pg/mL in plasma (citrate) and 896 pg/mL in plasma (heparin).

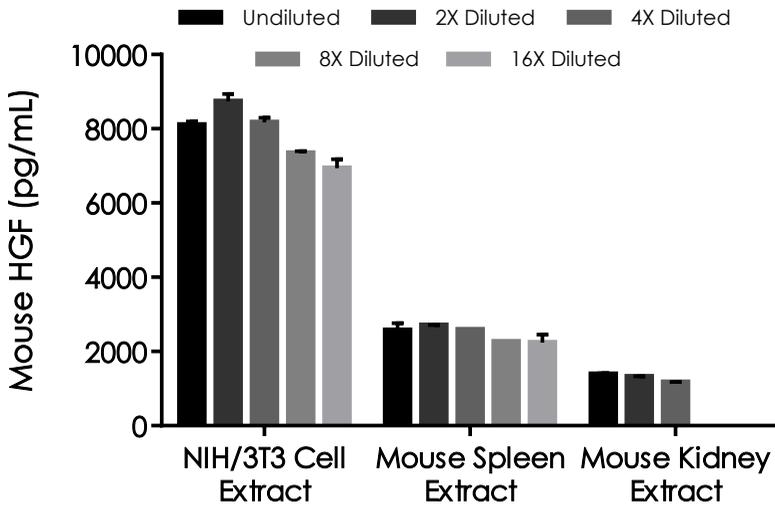

**Figure 4.** Interpolated concentrations of native HGF in Mouse cell and tissue extracts, extracts based on a 4,000  $\mu\text{g}/\text{mL}$  extract load. The concentrations of HGF were measured in duplicate and interpolated from the HGF standard curve and corrected for sample dilution. The interpolated dilution factor corrected values are plotted (mean  $\pm$  SD,  $n=2$ ). The mean HGF concentration was determined to be 7,859 pg/mL in NIH/3T3 cell extract, 2,472 pg/mL in mouse spleen tissue extract and 1,291 pg/mL in mouse kidney tissue extract.

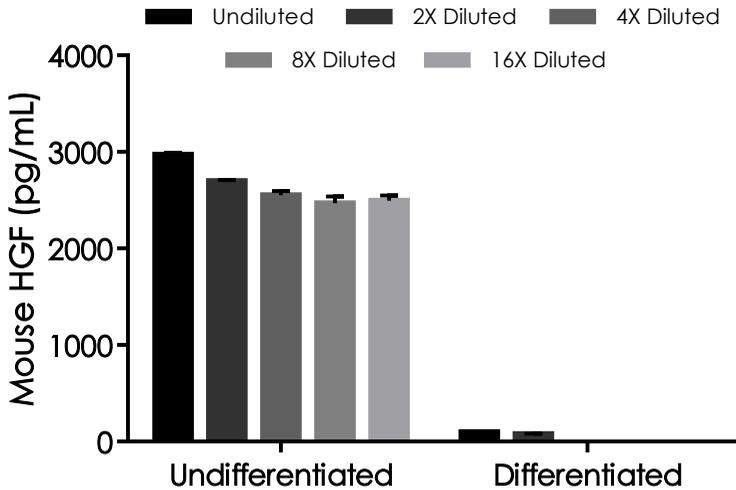

**Figure 5.** Interpolated concentrations of native HGF in Mouse 3T3-L1 undifferentiated and 3T3-L1 Adipocyte differentiated cell extract samples, sample based on a 1,000  $\mu\text{g/mL}$  extract load. Mouse 3T3-L1 cells were either undifferentiated (not allowed to reach confluency) or differentiated (cultured for 48 hours in adipocyte differentiation medium). The concentrations of HGF were measured in duplicate and interpolated from the HGF standard curve and corrected for sample dilution. The interpolated dilution factor corrected values are plotted (mean  $\pm$  SD,  $n=2$ ). The mean HGF concentration was determined to be 2,636 pg/mL in undifferentiated 3T3-L1 cell extract and 88 pg/mL in adipocyte differentiated 3T3-L1 cell extract.

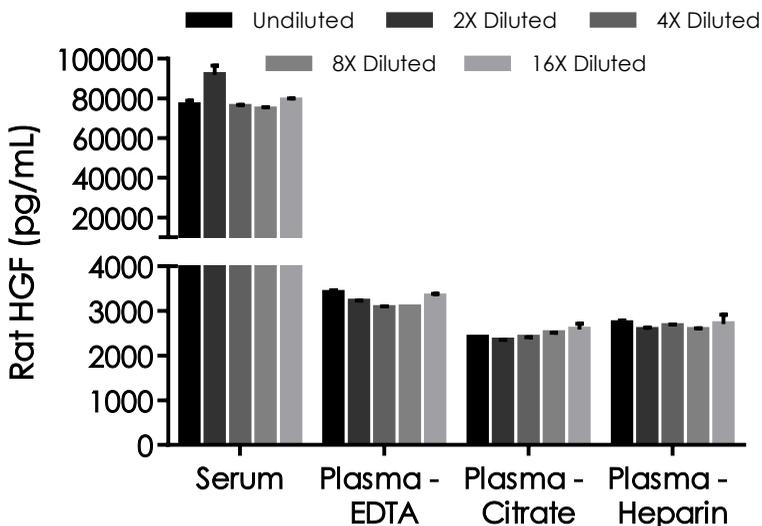

Fi

**Figure 6.** Interpolated concentrations of native HGF in Rat serum and plasma samples. The concentrations of HGF were measured in duplicates, interpolated from the HGF standard curves and corrected for sample dilution. Undiluted samples are as follows: serum 20%, plasma (EDTA) 100%, plasma (citrate) 100%, and plasma (heparin) 100%. The interpolated dilution factor corrected values are plotted (mean +/- SD, n=2). The mean HGF concentration was determined to be 79,543 pg/mL in serum, 3,151 pg/mL in plasma (EDTA), 2,395 pg/mL in plasma (citrate) and 2,611 pg/mL in plasma (heparin).

## 17. Assay Specificity

This kit recognizes both native and recombinant Mouse HGF protein in serum, plasma, and cell culture supernatant, cell and tissue extract samples only.

### **CROSS REACTIVITY**


Recombinant HGF Receptor was prepared at 4 ng/mL and assayed for cross reactivity. No cross-reactivity was observed.

### **INTERFERENCE**

Recombinant HGF Receptor was prepared at and 4 ng/mL and tested for interference. No interference with was observed.

## 18. Species Reactivity

This kit recognizes Mouse and Rat HGF protein.

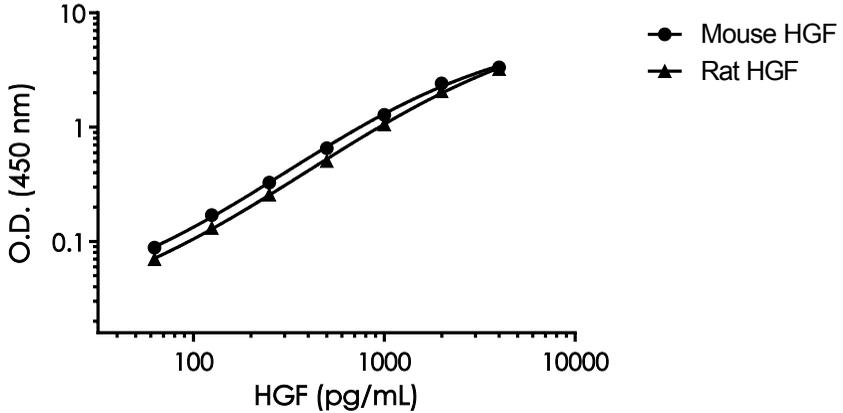

Recombinant Rat HGF was prepared at 4 ng/mL and diluted with a 2-fold titration alongside of Mouse HGF standard curve, and assayed for cross reactivity. 100% cross-reactivity was observed.

Recombinant Human HGF was prepared at 16 and 4 ng/mL and assayed for cross reactivity. No cross-reactivity was observed.

Please contact our Technical Support team for more information.

## 19. Troubleshooting

Problem	Reason	Solution
<b>Difficulty pipetting lysate; viscous lysate.</b>	Genomic DNA solubilized	Prepare 1X Cell Extraction Buffer PTR (without enhancer). Add enhancer to lysate after extraction.
<b>Poor standard curve</b>	Inaccurate Pipetting	Check pipettes
	Improper standard dilution	Prior to opening, briefly spin the stock standard tube and dissolve the powder thoroughly by gentle mixing
<b>Low Signal</b>	Incubation times too brief	Ensure sufficient incubation times; increase to 2 or 3 hour standard/sample incubation
	Inadequate reagent volumes or improper dilution	Check pipettes and ensure correct preparation
	Incubation times with TMB too brief	Ensure sufficient incubation time until blue color develops prior addition of Stop solution
<b>Large CV</b>	Plate is insufficiently washed	Review manual for proper wash technique. If using a plate washer, check all ports for obstructions.
	Contaminated wash buffer	Prepare fresh wash buffer
<b>Low sensitivity</b>	Improper storage of the ELISA kit	Store your reconstituted standards at -80°C, all other assay components 4°C. Keep TMB substrate solution protected from light.
<b>Precipitate in Diluent</b>	Precipitation and/or coagulation of components within the Diluent.	Precipitate can be removed by gently warming the Diluent to 37°C.

# 20. Notes





## Technical Support

Copyright © 2017 Abcam, All Rights Reserved. The Abcam logo is a registered trademark. All information / detail is correct at time of going to print.

### **Austria**

wissenschaftlicherdienst@abcam.com | 019-288-259

### **France**

supportscientifique@abcam.com | 01.46.94.62.96

### **Germany**

wissenschaftlicherdienst@abcam.com | 030-896-779-154

### **Spain**

soportecientifico@abcam.com | 91-114-65-60

### **Switzerland**

technical@abcam.com

Deutsch: 043-501-64-24 | Français: 061-500-05-30

### **UK, EU and ROW**

technical@abcam.com | +44(0)1223-696000

### **Canada**

ca.technical@abcam.com | 877-749-8807

### **US and Latin America**

us.technical@abcam.com | 888-772-2226

### **Asia Pacific**

hk.technical@abcam.com | (852) 2603-6823

### **China**

cn.technical@abcam.com | +86-21-5110-5938 | 400-628-6880

### **Japan**

technical@abcam.co.jp | +81-(0)3-6231-0940

### **Singapore**

sg.technical@abcam.com | 800 188-5244

### **Australia**

au.technical@abcam.com | +61-(0)3-8652-1450

### **New Zealand**

nz.technical@abcam.com | +64-(0)9-909-7829